

## Pipet-Lite XLS+ Single Channels

### Performance You Can Trust




#### Total hand comfort

With sure-fit handles, contoured finger hooks, lighter springs, low-drag seals and Rainin's patented Magnetic Assist™ technology, you can count on Pipet-Lite XLS+ to provide comfort all day long.



#### Smooth operation

The sealing system on XLS+ pipettes delivers a plunger stroke that is smooth, precise and requires less force than ever. Perfect for pipetting applications that require delicate control.

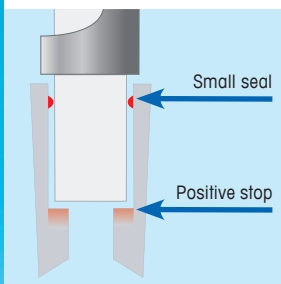

#### Lightest touch

Nothing beats Rainin's patented LTS LifeTouch™ System for ergonomics and ease of use. The small sealing area and positive stop ensure tips load quickly, seal perfectly and are easy to eject. And with new low-drag seals and lighter springs, your hands will feel the difference. No need to pound on the tips!



#### GLP/GMP security

Rainin pipettes were the world's first to be equipped with Radio Frequency Identification (RFID) chips for storing service and data profiles on-board, tracking calibration data and managing pipette inventory



### Pipet-Lite XLS+ with RFID

#### The Performance Champions

The new Pipet-Lite™ XLS+™ is one of the most ergonomic accurate and precise pipettes on the market. Its comfortable handle, light springs and "Magnetic Assist™" technology, ensure light and smooth operation, while significantly reducing the risk of repetitive strain injuries. Pipet-Lite XLS+ is the world's first pipette equipped with RFID tags for advanced calibration and asset management. Tip shaft options included low-force LTS™ for improved ergonomics and universal fit.

# Pipet-Lite™ XLS+™ with RFID

## The Performance Champion

### Feature Loaded

- Ergonomic handle design with comfortable fingerhook
- Improved volume adjustment lock prevents drift
- Built-in shock absorber inside the tip ejector reduces thumb impact
- Autoclavable shaft assembly
- Available with LTS and universal tip/shaft configurations
- Corrosion free, easy to clean tip ejector
- RFID tags for asset management



The optional RFID reader works with all Rainin RFID-enabled pipettes, and interfaces with METTLER TOLEDO's LabX™ Direct Pipette-Scan™ calibration check.

## Technical Data

### Specifications

Model	Volume (µL)	Increment (µL)	Accuracy		Precision	
			%	µL (±)	%	µL (±)
L-2XLS+, SL-2XLS+	0,2	0,002	12,0	0,024	6,0	0,012
	1		2,7	0,027	1,3	0,013
	2		1,5	0,030	0,7	0,014
L-10 XLS+, SL-10XLS+	1,0	0,02	2,5	0,025	1,2	0,012
	5,0		1,5	0,075	0,6	0,03
	10,0		1,0	0,1	0,4	0,04
L-20XLS+, SL-20 XLS+	2	0,02	7,5	0,15	2,0	0,04
	10		1,5	0,15	0,5	0,05
	20		1,0	0,2	0,3	0,06
L-100XLS+, SL-100XLS+	10	0,2	3,5	0,35	1,0	0,1
	50		0,8	0,4	0,24	0,12
	100		0,8	0,8	0,15	0,15
L-200XLS+, SL-200XLS+	20	0,2	2,5	0,5	1,0	0,2
	100		0,8	0,8	0,25	0,25
	200		0,8	1,6	0,15	0,3
L-300XLS+, SL-300XLS+	30	0,5	2,5	0,75	1,0	0,3
	150		0,8	1,2	0,25	0,375
	300		0,8	2,4	0,15	0,45
L-1000XLS+, SL-1000XLS+	100	2	3,0	3,0	0,6	0,6
	500		0,8	4,0	0,2	1,0
	1000		0,8	8,0	0,15	1,5
L-2000XLS+, SL-2000XLS+	200	2	3,0	6,0	0,6	1,2
	1000		0,8	8,0	0,2	2,0
	2000		0,8	16,0	0,12	2,4
*L-5000XLS, *SL-5000XLS	500	5	2,4	12,0	0,6	3,0
	2500		0,6	15,0	0,2	5,0
	5000		0,6	30,0	0,16	8,0
*L-10MLXLS, *SL-10MLXLS	1 mL	20	5,0	50,0	0,6	6,0
	5 mL		1,0	50,0	0,2	10,0
	10 mL		0,6	60,0	0,16	16,0
*L-20MLXLS	2 mL	20	5,0	100,0	0,6	12,0
	10 mL		1,0	100,0	0,2	20,0
	20 mL		0,6	120,0	0,16	32,0

### Ordering Information

Cat. No, LTS	MT-Order No, LTS	Cat. No, Traditional	MT-Order No, Traditional	Volume Range
L-2XLS+	17014393	SL-2XLS+	17014413	0,1–2 µL
L-10XLS+	17014388	SL-10XLS+	17014409	0,5–10 µL
L-20XLS+	17014392	SL-20XLS+	17014412	2–20 µL
L-100XLS+	17014384	SL-100XLS+	17014408	10–100 µL
L-200XLS+	17014391	SL-200XLS+	17014411	20–200 µL
L-300XLS+	17014405	SL-300XLS+	17014414	20–300 µL
L-1000XLS+	17014382	SL-1000XLS+	17014407	100–1000 µL
L-2000XLS+	17014390	SL-2000XLS+	17014410	200–2000 µL
*L-5000XLS	17011790	*SL-5000XLS	17011801	500–5000 µL
*L-10MLXLS	17011783	*SL-10MLXLS	17011795	1000 µL–10 mL
*L-20MLXLS	17011788	N/A	N/A	2000 µL–20 mL
L-STARTXLS+	17014406	SL-STARTXLS+	17014415	Starter Kits**

\* The 5000 µL, 10 mL and 20 mL models do not appear the same as the photos.

\*\* Include 20, 200, 1000 µL models, tips and accessories.

### Accessories

Cat. No,	MT-Order No,	Description
RFID-KIT	17011966	RFID Kit including Rainin RFID Reader and LabX Direct Pipette-Scan Software
RFID-RDR	17011964	Rainin RFID Reader
LABX-PIPET	17011965	LabX Direct Pipette-Scan Software

Rainin, Pipet-Lite, XLS, XLS+, LTS, Pipetting 360°, LiteTouch, and Magnetic Assist are trademarks of Rainin Instrument, LLC, Direct Pipette-Scan and LabX are trademarks of METTLER TOLEDO Corporation.



### Rainin Instrument, LLC

7500 Edgewater Drive, Oakland, CA 94621  
Phone +1 510 564 1600

a METTLER TOLEDO Company

Subject to technical changes  
© 04/2014 Rainin Instrument, LLC  
Printed on demand from online content. 17700502

[www.mt.com/rainin](http://www.mt.com/rainin)

For more information